


## Foley Catheter

Product Model : **Default Model No**

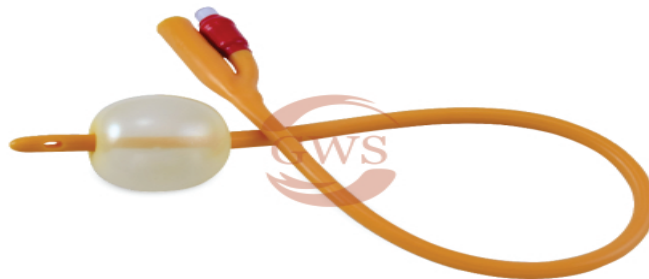


**Category: Medical Equipment, Medical Disposables, Urology,**

### Summary :

Our '**Foley Catheter**' is a flexible tube made from natural latex rubber, that a clinician passes through the urethra and into the bladder to drain urine. It is the most common type of indwelling urinary catheter. The tube has two separated channels, or lumens, running down its length. One lumen, open at both ends, drains urine into a collection bag. The other has a valve on the outside end and connects to a balloon at the inside tip. The balloon is inflated with sterile water when it lies inside the bladder to stop it from slipping out. The catheter is a 16F unit which measures at 5.3mm in diameter.

With ISO 9001-2015; ISO 13485-2016 certification and WHO-GMP compliance we can ensure the best quality and medical safety of our product.

**Legacy of Excellence in Healthcare...**

 [info@gwsmed.com](mailto:info@gwsmed.com), [gwssurgicals@gmail.com](mailto:gwssurgicals@gmail.com)

 +91 124 4915548  [www.gwsmed.com](http://www.gwsmed.com)



**Corp Office:** 1205, 12<sup>th</sup> Floor Ocus Quantum Sector 51, Gurugram - 122003 (Hr.) India

**Works:** Plot No. 38, Sector 7, Sidcul Industrial Area, Haridwar-249403 (Uttarakhand), India

