


Five Functional ICU Bed Electric (ABS Panel & ABS Side Railing)

Product Model : **HF-101**

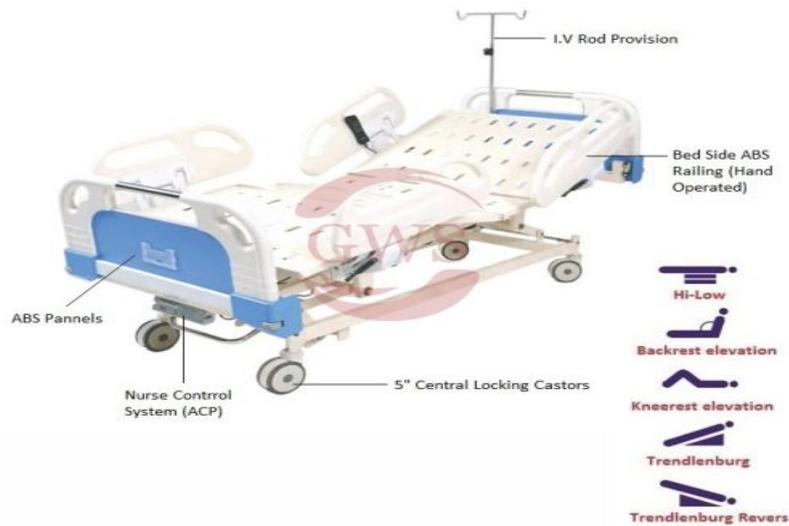

Category: hospital furniture, Medical Equipment, ICU Beds,

Legacy of Excellence in Healthcare...

 info@gwsmed.com, gwssurgicals@gmail.com

 +91 124 4915548  www.gwsmed.com



Corp Office: 1205, 12th Floor Ocus Quantum Sector 51, Gurugram - 122003 (Hr.) India

Works: Plot No. 38, Sector 7, Sidcul Industrial Area, Haridwar-249403 (Uttarakhand), India

Description

- Advance Linear Actuator System [01pc. Motor + 01pc. Control Box + 01pc. Nurse Control System (ACP)]
- Four Individual Side Railing, ABS Moulded.
- ABS Head and Foot Panel with Locking and Corner Buffer.
- M.S Tubular Framework.
- 4 Section Top made-up of Oval perforated CRCA M.S Sheet.
- Bed Mounted on 5" Central Locking Castors.
- S.S Telescopic IV Rod Provisions on 4 Locations.
- Urine bag holder.
- Sitting Position Adjustable.
- Pre-treated & Epoxy Powder Coated.
- Overall Size (apprx.): 2090L x 910W x 460-700H mm
- Functions by Hand Remote: (a) Hi-low (b)Backrest elevation (c)Kneerest elevation (d)Trendlenburg / Reverse Trendlenburg
- Functions by Nurse Control System (ACP): (a)Hi-low with locking function for handset. (b)Backrest elevation with locking function for handset. (c)Kneerest elevation with locking function for handset. (d)Trendlenburg / Reverse Trendlenburg with locking function. (e)Electric CRP for Zero Positioning of Bed.