


Portable Anaesthesia Machine

Product Model : **AN-164**

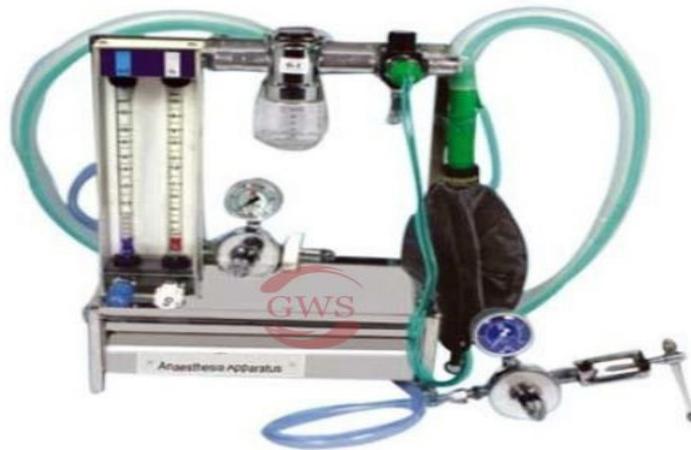

Category: Anaesthesia Products & Equipment, Anaesthesia Machine,

Legacy of Excellence in Healthcare...

 info@gwsmed.com, gwssurgicals@gmail.com

 +91 124 4915548  www.gwsmed.com



Corp Office: 1205, 12th Floor Ocus Quantum Sector 51, Gurugram - 122003 (Hr.) India

Works: Plot No. 38, Sector 7, Sidcul Industrial Area, Haridwar-249403 (Uttarakhand), India

Description

- Construction: trolley made of tubular rigid steel section electrostatically powder coated or stainless steel.
- Suitable drawer
- Unitised rotamotor block for oxygen & nitrous oxide with fine adjustment controls & colour coated knob and luminous back plate.
- Oxygen flow meter calibrated from 100cc/mm to 100 ltrs. / Minute & Nitrous oxide flow meter from 200cc/min to 12 liters/min down Stream receipt of oxygen
- One gold man type vaporizer for halothene
- Non return cum emergency oxygen unit with blow off valve in circuit ensures that pressure does not build-up beyond 200cm. water column & minimises risk of back flow of gases.
- Two pre set regulators (forged quality) one each for O₂ & N₂o with high pressure tubing non interchangeable for connection with flow meter inlet
- One Maggils & brain circuit